# concept

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## ACKNOWLEDGEMENT

Thank you for purchasing a Concept product. We hope you will be satisfied with our product throughout its service life.

Please study the entire Operating Manual carefully before you start using the product. Keep the manual in a safe place for future reference. Make sure other people using the product are familiar with these instructions.

TECHNICAL PARAMETERS	
Voltage	230 V ~ 50 Hz
Power input	12 W

## IMPORTANT SAFETY PRECAUTIONS:

- Make sure the mains voltage corresponds to the values stated on the product's rating plate.
- The unit may not be cleaned, stored, or covered until it has cooled down completely.
- Do not leave the unit unattended while it is turned on or even only plugged into the wall outlet.
- When plugging the unit into the wall outlet and disconnecting it from power network the mode selector has to be switched to **O** (off).
- When disconnecting the unit from the wall outlet, never pull the power cord. Grasp the plug and disconnect it by pulling.
- Do not allow children or unskilled people to handle the unit. Use the unit out of the reach of these individuals.
- Take extra care when using the unit **near children**.
- Do not cover the unit, or use it under a sheet or duvet – there is a risk of overheating.
- Use the unit only for the intended purpose specified in the Operating Manual.
- Make sure the power cord does not touch hot surfaces.
- **Turn the unit off, unplug the power cord from the wall outlet, and let the unit cool down before cleaning, assembly and disassembly, and immediately after use.**
- Do not use abrasive or chemically aggressive substances to clean the unit.
- Do not use the unit if it does not operate properly, if it has been dropped, damaged or exposed to any liquids. Refer the unit to an authorized service center for testing and repair.
- The unit **is not designed for outdoor use**.
- The unit is suitable **for domestic use only. It is not designed for commercial operation**.
- Use only accessories recommended by the manufacturer.
- Do not touch the unit with wet or damp hands.
- **Do not immerse the power cord, plug, or unit in water or any other liquid.**
- Never use an appliance with a damaged power cord or plug. Contact an authorized service center immediately to have the defective components repaired or replaced.
- Do not repair the appliance yourself. Contact an authorized service center.
- This appliance has a heated surface. People who do not feel heat have to be careful when using the appliance.
- Do not use the unit near a shower, bathtub, washbowl, swimming pool, or other vessels containing water.
- Do not use the unit near the eyes, and do not look directly in the infrared beams.


Meaning of the label: The unit is not designed for use in the bathtub or shower.

**Failure to follow the manufacturer's instructions may cause refusal of warranty repair.**

## DESCRIPTION OF THE PRODUCT

1. **Massage head with infrared lamp**
2. **Massage surface**
3. **Flexible turning neck**
4. **Joint bending button**
5. **Mode selector**
- 6, 7. **Body, neck and head massage extensions**
8. **Muscle massage extension**
9. **Massaging brush**
10. **Supply cable**
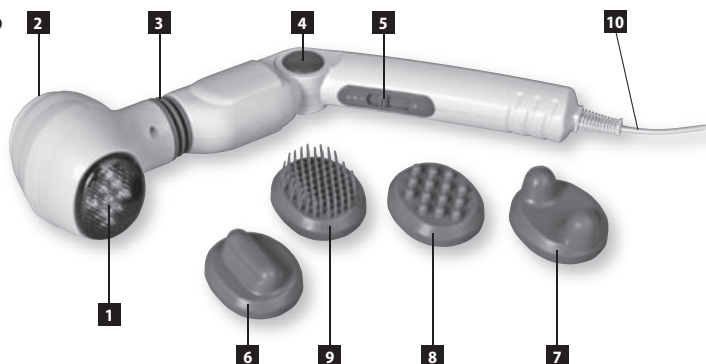
### NOTE

- Pregnant women, diabetics, people suffering from varicose veins and ill people should consult their doctor prior to use.
- A massage should be a pleasant and relaxing experience. If you feel pain, or any discomfort, consult your doctor.
- Do not use the unit on irritated, burnt or swollen places, injured skin or a rash. Avoid direct effects of the massage on joints, spine, front of the neck and face. If you have a cardiac stimulator installed, avoid the heart area.
- Do not use the unit for more than 20 minutes. Let it cool down for at least 20 minutes after use.

Prior to using a new unit, wipe the surface with a wet cloth and wash the accessories in warm water for hygienic reasons.

## OPERATING INSTRUCTIONS

1. Uncoil the supply cable.
2. Connect the plug of the power cord to the wall outlet.
3. You may adjust the shape of the unit according to massaged spot. The massaging head is rotary, and the handle's angle can be adjusted by pressing the joint button.
4. Select the required massage, or infrared mode.  
Mode selector positions:
  - 0 – Off.
  - 1 – Massage – Use the extension on the massaged spot according to the required massage mode.
  - 2 – Infrared lamp – prior to use on the selected spot, check the massage head temperature. The temperature can be regulated by switching the unit on and off.
  - 3 – Massage and infrared lamp – prior to use on the selected spot, check the massage head temperature.
5. Do not use the unit for more than 20 minutes. Let it cool down for at least 20 minutes after use.
6. Switch the mode selector to the **0** position (off), and disconnect the unit from the wall outlet after use.

### Information on the use of the infrared lamp

1. Therapeutic feature: The heat of the infrared beams reduces the pain of rheumatic and muscle diseases. Consult with your doctor prior to use.
2. Cosmetic feature: The infrared radiation supports blood circulation, and helps to open the pores. The skin absorbs cream and other cosmetic preparations better after application. Put the massage head with infrared lamp lightly on the skin, and move it gently so that the beams do not point to one spot only. You may repeat the application 3-4 times a day for 5-10 minutes, with at least 45-minute breaks. It is common for the skin to become slightly red during the application. If you feel discomfort, or the skin becomes too red, do not continue the application.

### Information on the massage feature

Move the massaging head in light circular motions on the selected spot for 6 to 8 minutes. Avoid the continuous massage of just one spot. The unit can be used through clothing, but direct application on the body ensures better results.

By selecting a proper extension you may achieve higher intensity of the massage.

Body, neck and head massage extensions (6, 7) will provide intensive massage of sore muscles, and at the same time they are convenient for massaging the feet. Muscle massage extension (8) helps to reduce and attenuate pain and tension in the muscles. The massage brush (9) performs intensive skin massage, and simultaneously supports blood circulation. It may be used to massage the scalp to promote better blood circulation to the hair roots.

## MAINTENANCE AND CLEANING

### WARNING!

Turn off the unit, unplug the power cord from the wall outlet, and let the unit cool down before cleaning, assembly and disassembly, and immediately after use. Prior to handling, make sure the unit has cooled down!

The surface of the unit may only be cleaned using a wet cloth. Do not use any detergents or hard objects as they may damage the unit's surface!

**Never clean the unit under running water, do not rinse it, and do not submerge it in water!**

## REPAIRS & MAINTENANCE

Any extensive maintenance or repair requiring access to the internal parts of the product shall be performed by an authorized service center.

### ENVIRONMENTAL PROTECTION

- Packaging materials and obsolete appliances should be recycled.
- The transport box may be disposed of as sorted waste.
- Polyethylene (PE) bags shall be handed over for recycling.

#### Recycling at the end of the service life:

A symbol on the product or its packaging indicates that this product should not be disposed of as household waste. It must be taken to the collection point authorized for recycling electric and electronic equipment. By making sure this product is disposed of properly, you will help prevent negative effects on the environment and human health that would otherwise result from inappropriate disposal of this product. You can learn more about recycling this product from your local authorities, a household waste disposal service, or in the shop where you purchased this product.

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