concept

Sonický elektrický zubní kartáček dětský
Sonická elektrická zubná kefka pre deti
Soniczna szczoteczka elektryczna
dla dzieci
Szónikus elektromos fogkefe gyerekek-
nek
Sonic elektriskā zobu birste bērniem
Sonic electric toothbrush for children

Elektrische Schallzahnbürste für Kinder
Brosse à dents électrique sonore pour
enfants
Spazzolino elettrico sonico per bambini
Cepillo eléctrico sónico para niños
Periuță de dinți electrică sonică pentru
copii

ZK3010 / ZK3020

CZ

SK

PL

HU

LV

EN

DE

FR

IT

ES

RO

Poděkování

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Acknowledgements

Thank you for purchasing a Concept product. We wish you much joy with your new appliance every day you use it.

Please read the instruction manual carefully before initial use. To refer to this manual any time you need to, we recommend you to keep it in a safe place. And please pass it to any future owner of the appliance.

TECHNICAL SPECIFICATIONS	
Voltage	5 V
Input	1,2 W
Function	Clean, Soft
Oscillation	24 000 strokes/min
Charging	Magnetic USB charger
Charging time	4 hours
Battery	DC3.7V ,260mAh lithium battery
Water resistance	IPX7
Operating temperature	5-40 °C

IMPORTANT SAFETY PRECAUTIONS

- Do not use the device in any other way than as described in this manual.
- If you have had dental surgery in the last two months, consult your dentist about using the device.
- During the first days of using the device, you may notice a slight bleeding from the gums. If the gums bleed for more than 14 days, stop using the device and consult your dentist.
- Do not use brush heads other than those recommended by the manufacturer.
- If you have a pacemaker or other device installed, consult your doctor about using this toothbrush.
- The device is intended only for cleaning teeth, gums and tongue. Do not use it for other purposes.
- If you experience discomfort or pain, stop using the device and seek medical attention.
- Remove all the covering and marketing materials from the device before the first use.
- Make sure the mains voltage corresponds to the values on the rating plate of the device.

- This product is not a toy. Children should be supervised at all times when the appliance is switched on to ensure that they are using the toothbrush correctly and that they do not play with the toothbrush.
- If you put the device down during use, always switch it off.
- Before connecting or disconnecting the device from the mains socket, make sure that the device is switched off.
- When disconnecting the device from the wall outlet, never pull the power cord, grasp the plug end and disconnect it by removing it from the receptacle.
- Do not use the device in the bathtub.
- Do not leave the device in a place from which it can fall into a sink or a bathtub.
- Do not soak the product in water for more than 30 minutes.
- Do not share the brush head with other people.
- Do not place metal objects (coins, staples, etc.) on the charger.
- Persons with reduced mobility, reduced sensory perception, mental disadvantage or persons not familiar with the operation must use the device only under the supervision of responsible person, who is familiar with the instructions.
- Pay extra attention if the device is used near children.
- Do not use the device as a toy.
- **Before cleaning and after use of the device, switch the device off and disconnect it from the mains.**
- Do not let the power cord touch hot surfaces.
- Do not use the device in an environment where there are explosive gases and flammable substances (solvents, varnishes, adhesives etc.).
- Do not use coarse and chemically aggressive substances to clean the device.
- Do not use the device if it is not working properly, has been dropped, damaged or immersed in liquid. Have it checked and repaired by an authorized service centre.
- Do not use the device with a damaged cord or plug and have it repaired immediately by an authorized service centre.
- Do not use the device outdoors.

- Avoid placing in high temperature conditions or in direct sunlight.
- The device is only suitable for home use and is not intended for commercial use.
- **Do not immerse the power cord, plug or device in water or any other liquid.**
- Do not use any accessories other than those recommended by the manufacturer.
- Do not disassemble, repair or modify the unit or parts of this equipment yourself. Contact an authorized service centre.
- This device can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the device in a safe way and understand the hazards involved. Cleaning and maintenance performed by the user must not be carried out by children unless they are 8 years old and under supervision. Children under the age of 8 must be kept away from the device and its power supply. Do not let children play with the device.

Damages from failure to follow the manufacturer's instructions are not covered by the warranty.

CAUTION:


This device is waterproof. It is suitable for cleaning under running water. Do not use or clean the device while it is connected to the power cord.

PRODUCT DESCRIPTION

- 1 Replaceable brush head
- 2 Led Observation Light/Colour Light
- 3 On/Off button and mode selection
- 4 Handle
- 5 Magnetic USB charging port

Package contains:

- 1x Rechargeable toothbrush
- 2x Smaller brush head (1-2 years)
- 2x Larger brush head (3-5 years)
- 1x Charging cable
- 1x User manual"

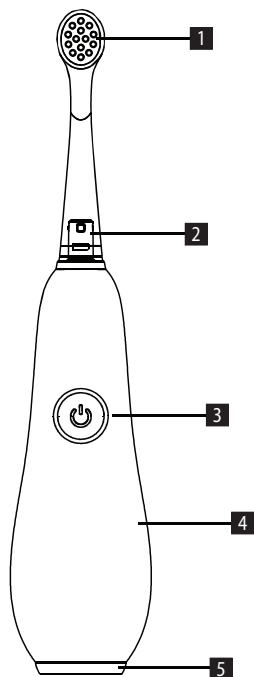OPERATING INSTRUCTIONS

Charging


- Fully charge the battery before using the device for the first time.
- If the device has not been used for more than 6 months, charge it for more than 15 hours.
- If the battery capacity is less than 20%, the white LED will flash.
- Charge only when the battery is low.
- To make better use of the battery, we recommend discharging the battery once a year.
- Connect the device to the magnetic charging port. Plug the USB cable into the power adapter. The LED light will flash while charging (full charge will occur in about 4 hours). When fully charged, the charger will automatically shut off (overcharge protection) and the LED light will go out.
- If you are using the device for the first time or if you have not used it for a long time: the indicator may not light up during the first 5 min of charging. This means that the battery is almost empty, which is a normal phenomenon.

Smart timer


- The recommended cleaning time is at least 2 minutes. The device switches off automatically after two minutes.
- Every 30 seconds, the device pauses to alert you to move to the next part of your mouth.

Cleaning modes


- This device has 2 modes: **Clean**, **Soft**. It helps to remove plaque more effectively, to have healthier teeth and gums.
- **Clean** - standard choice for cleaning, removes 100 % of plaque
- **Soft** - gentle cleaning of sensitive teeth and gums

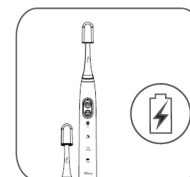
Fully charge the battery before using the device for the first time - LED light goes out.

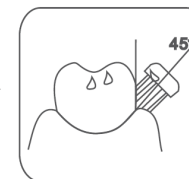
To clean your gums effectively, perform the cleaning at an angle of 45°.

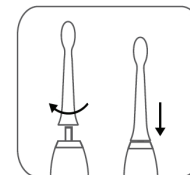
Place the brush head on the brush body. To remove the brush head, reverse the procedure.

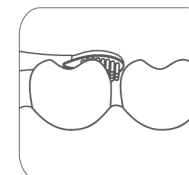Apply the brush lightly to avoid damaging the gums.

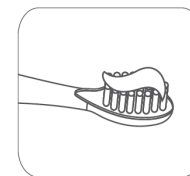Rinse the brush head, apply the paste, put the brush head in your mouth and then turn on the brush.

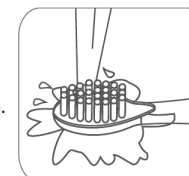Rinse the brush head with water and dry.

TROUBLE SHOOTING

Problem	Possible cause	Remedy
The strength of the vibration fades	Excessive pressure exerted on your teeth.	Carefully place the brush head on the tooth.
Loose brush head or abnormal sound	The brush head is damaged.	Replace the brush head.
The running time of the brush is short	The handle is not positioned correctly, leading to charging failure.	Make sure the handle is correctly positioned.
	Low battery.	Extend the charging time.
The toothbrush does not work	Low battery capacity. Other cause.	The battery needs to be replaced. Please contact customer service to resolve the problem.

CLEANING AND MAINTENANCE

Brush head - After use, rinse with water and dry.

Brush body - Rinse with water and dry. Do not use water warmer than 40 °C or any cleaning agents.

Charging base - Wipe only with a dry cloth. Wet and damp fabrics can cause a short circuit.

Replacing the brush head - Dentists recommend replacing the brush head every 3 months.

Store the device in a dry and dust-free place.

Do not expose to high temperatures.

Keep out of reach of children.

Warning: The device contains a rechargeable lithium battery. Recycle accordingly.

SERVICE CENTRE

Extensive maintenance or repair, which requires access to internal parts of the product, has to be performed by a qualified service centre.

ENVIRONMENTAL CONCERNS

- Packaging materials and old appliances should be recycled.
- Packaging materials may be disposed of as sorted waste.
- Dispose the plastic bags made of polyethylene (PE) of as sorted waste.

Recycling of the appliance at the end of its lifespan

This appliance is labelled in accordance with European Directive 2012/19/EU concerning used electrical and electronic appliances (waste electrical and electronic equipment - WEEE). The symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The disposal of such material has to be done according to the recycling regulations. For more detailed information about recycling of this appliance see the authorized local office or service for household waste disposal or the shop, where you purchased the appliance.

The product meets all the necessary requirements of the EU directives applicable to the specific product.

Changes in text, design, and technical specifications are subject to change without notice and we reserve the right to change them.

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